OMRON

Manual Inflation Blood Pressure Monitor

MODEL HEM-4030

Slim, Compact design for easy measurement and storage



- Hypertension Indicator (when systolic/ diastolic is out of standard range of 135/85mm Hg)
- Memory - 21 sets of readings
- Durable, easy to use inflation bulb


eco
No Pollutants

A Good Sense of Health

TRUSTED WORLDWIDE

Omron HEM-4030 Manual Inflation Blood Pressure Monitor uses the oscillometric method of blood pressure measurement. It comes with an easy to use inflation bulb for quick and easy manual inflation.

Hypertension Indicator

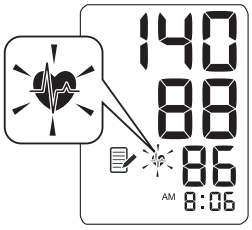


If your systolic or diastolic presure is outside the standard range, the heartbeat symbol will blink when the measurement result is displayed. Recent research suggests that the following values can be used as a guide to high blood pressure for measurments taken at home.


Systolic Blood Pressure	Above 135 mmHg
Diastolic Blood Pressure	Above 85 mmHg

This criteria is for home blood pressure measurement.

Correct usage and cuff application

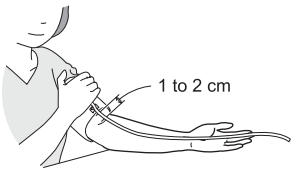

- The bottom edge of the cuff should be 1 to 2 cm above the elbow.
- Marker (arrow under tube) is centred on the middle of your inner arm.

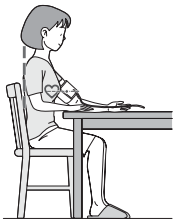

- Sit on a chair with your feet flat on the floor.
- Sit upright with your back straight.
- The cuff should be at the same level as your heart.

Specifications:

Model	HEM-4030
Display	LCD Digital Display
Measurement Range	Pressure: 0-299 mmHg Pulse: 40-180/minute
Accuracy	Pressure : ±3 mmHg Pulse : ±5% of reading
Measurement method	Oscillometric method
Inflation	Manual by Inflation bulb
Deflation	Automatic pressure release valve
Pressure Detection	Capacitive Pressure Sensor
Power source	2 “AAA” batteries 1.5V
Battery life	Approx. 1500 uses when used twice a day
Operating Temperature	10°C to 40°C, 30% to 85% RH
Storage Temperature	-20°C to 60°C, 10% to 95% RH
Console weight	Approx. 75g (excluding batteries)
Unit Dimensions	Approx. 105 (L) x 64 (W) x 35 (H) mm
Cuff Dimensions	Approx. 146 (W) x 446 (L) mm (Medium cuff: fits arm circumference 22-32cm)
Cuff Material	Nylon and polyester
Contents	Main unit, Cuff, Inflation bulb, Instruction Manual, storage case, batteries
Optional	S size fan cuff (fits arm circumference 17-22cm), M size fan cuff (fits arm circumference 22-32cm), L size fan cuff (fits arm circumference 32-42cm), Small cuff and bulb combination, Regular bulb

Note: Should a small cuff is required, please note that the small cuff and small bulb should be used together.
Specifications are subject to change without prior notice.

OMRON HEALTHCARE Co.,Ltd
24 Yamanouchi Yamanoshita-cho.
Ukyo-ku, Kyoto 615-0084 JAPAN

Distributed by:

Asia Pacific HQ
OMRON HEALTHCARE SINGAPORE PTE LTD
438A Alexandra Road, #05-05/08,
Adexandra Technopark, Singapore 119967
www.omron-healthcare.com.sg

EC REP
OMRON HEALTHCARE EUROPE B.V.
Kruisweg 577, 2132 NA Hoofddorp
THE NETHERLANDS / Les Pays-Bas

CE 0197