
beurer

HK 67

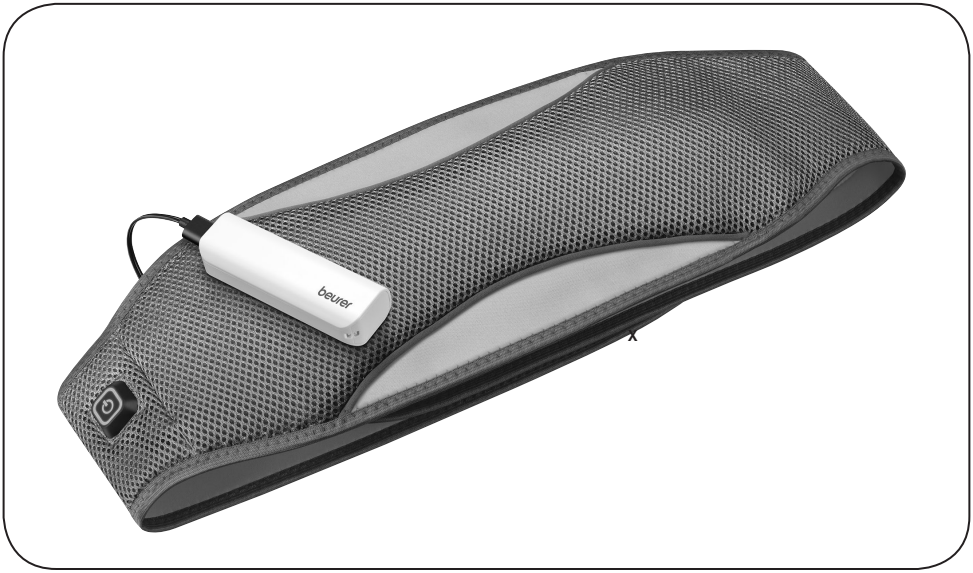READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT
Keep this manual in a safe location for future reference

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ENGLISH **Mobile Heating Belt To Go**
Instructions for use 2

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ENGLISH

CONTENT

1. IMPORTANT SAFETY INSTRUCTIONS 2
2. Package Contents 4
3. Parts and Controls 4
4. Using the Power Bank 4
5. To Use 5
6. Cleaning and Maintenance 6
7. Troubleshooting Guide 6
8. Technical data 7
9. FCC Compliance Statement 7
10. Warranty 7

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

Dear Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Please read these instructions for use carefully and keep them for later use. Make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer team

1. IMPORTANT SAFETY INSTRUCTIONS

Signs and symbols

The following signs appear in the Safety Section (page 2) and in this manual on page 3, 4, 5, 6.

READ AND UNDERSTAND THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

	Read all instructions before use		Do not bleach
	Keep away from sharp objects		Do not tumble dry
	Do not use when folded or bunched up		Do not iron
	Keep away from children under three years of age		Do not dry clean
	Machine wash warm, gentle cycle		

⚠ **WARNING**	**WARNING:** Indicates a hazardous situation which, if not avoided, could result in death or serious injury.
⚠ **CAUTION**	**CAUTION**: Indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.
ⓘ **NOTICE**	**NOTICE:** addresses practices not related to personal injury, such as product and/or property damage.

⚠ WARNING – To reduce the risk of fire, electric shock, or serious personal injury:

- This heating belt is only designed to warm up the human body. Do not use this product for any other purpose.
- Non-observance of the following notes may result in personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information is not only intended to protect your health and the health of others, it should also protect the product. For this reason, pay attention to these safety notes and include these instructions when handing over the product to others.
- This heating belt must not be used by persons who are not sensitive to heat or by vulnerable persons who may not be able to react to overheating.
- This heating belt must not be used by children under three.
- This heating belt is intended for adult use and should not be used by children unless the controller has been set by a parent or responsible adult and the child is closely supervised at all times during use.
- This heating belt may not be used by persons (including children) with reduced physical, sensory or mental skills or a lack of experience or knowledge, unless they are closely supervised and have been instructed on how to use the heating belt safely and are fully aware of the risks of use.
- This product is not a toy. Children must not play with it.
- Cleaning and user maintenance must not be performed by children.
- This heating belt is not designed for use in hospitals.
- This heating belt is suitable for personal use only, not for medical or commercial purposes.
- Do not insert sharp objects into the product.
- Do not use when folded or bunched up.
- Do not use the product if you or the product is wet.
- This heating belt may only be used in conjunction with the Power Bank type WE-213 provided with the unit.
- Do not expose the Power Bank to direct sunlight or high temperatures.
- The electrical and magnetic fields emitted by this heating belt may interfere with the function of a pacemaker. However, they are still well below the limits for electrical field strength (5000 V/m max), magnetic field strength (80 A/m max), and magnetic flux density (0.1 millitesla max). Please consult your doctor and the manufacturer of your pacemaker before using this heating belt.
- Do not pull, twist or make sharp bends in the cables.
- Please check this heating belt and connection cable frequently for signs of wear and tear or damage. If detected, discontinue use and contact customer service.
- When this heating belt is switched on:
 – Do not place any objects on it (e.g. suitcase or washing basket).
 – Heat sources, such as hot water bottles, heating belts or similar, must not be placed on the belt.
- While charging, do not cover the Power Bank with covers, cushions or similar items.
- Do not charge the Power Bank while the heating belt is in use.
- Do not charge the Power Bank using a computer as it could become overloaded.
- Risk of explosion! Do not put the Power Bank into fire or it may explode.
- Do not disassemble or crush the Power Bank.
- Keep the Power Bank away from direct sunlight and high temperatures.
- Always observe the notes on operation (chapter 4), cleaning, maintenance and storage (chapter 5).

- Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding correct charging must be observed at all times.
- For pluggable equipment, the socket-outlet shall be installed near the equipment and shall be easily accessible.
- Never expose the power bank to extreme heat. This applies in particular to storage in a car. Extreme temperatures can occur inside the vehicle and in the glove compartment when the vehicle is parked for a long period of time. Therefore, remove the power bank from the vehicle.
- If you have any questions about using this device, contact Customer Service.


⚠ **CAUTION** This heating belt is equipped with a Safety Shut-off Sensor that prevents the heating belt from overheating. If the Safety Shut-off Sensor is triggered, the Power Button on the heating belt and the white LED on the power bank remain illuminated. After cooling down, the heating belt switches on automatically.

2. Package Contents

1 x Mobile Heating Belt To Go HK67
1 x Power Bank Type WE-213
1 x Charging Cable
1 x Storage Bag
1 x Instruction Manual

3. Parts and Controls

1. Power Bank Type WE-213
2. Heating Belt
3. Power Button
4. Power Cable
5. Charging Cable
6. Storage Bag

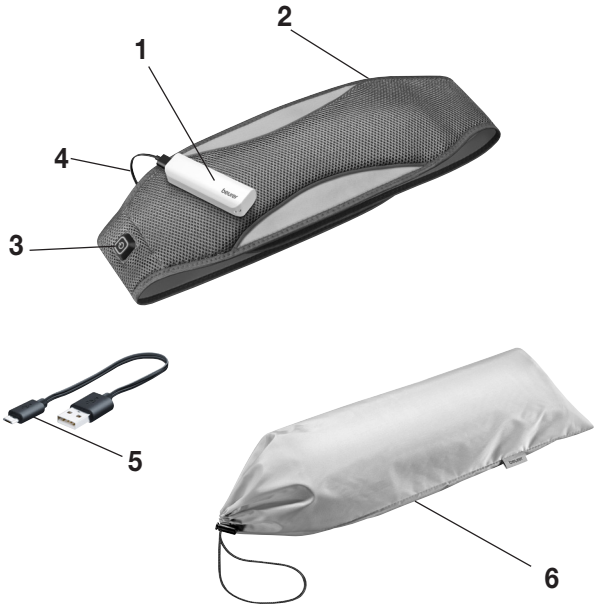4. Using the Power Bank

Charge the Power Bank before initial use and if the device has not been used for a long period.
To charge:

1. Insert the provided Charging Cable into the connection socket of the Power Bank.
2. Connect the other end of the Charging Cable to a USB charger via the USB port (USB charger not included) with at least a 1-Amp output. The Battery Level LED on the Power Bank will flash white while charging. The charging/decharging levels are as follows:
 1 x flash: 0 % - 25 %
 2 x flash: 26 % - 50 %
 3 x flash: 51 % - 75 %
 4 x flash: 76 % - 99 %
 Steady on 100 %
3. The Power Bank is fully charged when the Battery Level LED stops flashing and displays steadily (about 4.5 hours for a completely discharged Power Bank).

4. Disconnect the Charging Cable; the LED will shut off. Charge the Power Bank at least once per year to maintain its capacity.

WARNING

- While charging, do not cover the Power Bank. Only charge the Power Bank in a dry environment at normal room temperature.

Charging other devices from the Power Bank

The Power Bank can also be used to charge other rechargeable devices, such as smartphones. To charge your device using the Power Bank, simply connect your device's USB cable to the Power Bank and to the USB port on your device.

5. To Use

1. Connect the integrated Power Cable on the heating belt with the USB Power Socket on the Power Bank; the Power Button will flash white three times, indicating the device is ready for use.
2. Place the heating belt around your waist and secure the ends together using the hook and loop fasteners.

CAUTION

- Ensure that the heating belt does not bunch up or fold over during use.

3. Press the Power Button on the heating belt to activate; the Power Button will illuminate in green indicating the device is on Low. The power bank's LED illuminates or flashes white, according to the charging level of the power bank.
4. To increase the heat level, press the Power Button again; the Power Button will illuminate yellow to indicate the medium setting. Press a third time to increase to maximum heat; the Power Button will illuminate red.
5. To shut off the heating belt, press the Power Button again. The LED on the heating belt will shut off immediately, the LED on the power bank after a few seconds.

WARNING

- If the heating belt is being used over several hours, set the lowest temperature setting on the Power Bank to avoid overheating the heated body part, which may lead to burns to the skin.

NOTICE: Protect the Power Bank from extreme cold to maintain its performance.

NOTE

- When the Power Bank level is too low, the heating belt will shut off. The Power Bank must then be charged before using it again.

Setting the temperature

Setting 1: Minimum heat (LED display of the On/Off button lights up green)
Setting 2: Medium heat (LED display of the On/Off button lights up yellow)
Setting 3: Maximum heat (LED display of the On/Off button lights up red)
OFF: On/Off button not lit

Automatic shut-off
This heating belt is equipped with an automatic shut-off function that shuts the heating belt off after about 110 minutes. The Power Button of the heating belt and the white LED on the Power Bank will shut off. Press the Power Button to switch the heating belt back on. If the heating belt does not switch on, the Power Bank is depleted and must be recharged.

6. Cleaning and Maintenance

⚠ WARNING

- Prior to cleaning, always disconnect the Power Bank from the heating belt and remove it from the pouch.
- Always keep the Power Bank away from water or other liquids, as this may cause damage. To clean the Power Bank, use a dry, lint-free cloth. Do not use any chemical or abrasive cleaning agents.

NOTICE: To prevent product damage, do not dry clean, wring out, tumble dry, or iron the heating belt. To machine wash, place it in a pillow case. Set the washing machine to the delicate or wool cycle and the warm water setting. Always wash the heating belt together with other items. Use a mild laundry detergent according to the manufacturer's instructions.

Frequent washing of the heating belt is not recommended and should only be machine washed a maximum of 10 times during its life. Let the heating belt dry flat on a drying rack. To prevent product damage, do not use pegs or similar items to attach the heating belt to the drying rack.

NOTICE: Reconnect the heating belt to the Power Bank only after the heating belt and the connector plug are completely dry.

Storage

A storage bag is provided for taking the heating belt with you on journeys. Remove the Power Bank. Carefully roll up the heating belt and then also place the Power Bank in the storage bag. If you do not plan to use the heating belt for a long period of time, store it in the original packaging in a cool, dry place. Always remove the Power Bank before storing the device.

The Power Bank should be charged at least once per year to maintain its capacity.

NOTICE: Allow the heating belt to cool down before storing. To avoid making sharp bends in the heating belt, do not place any objects on top of it while storing.

7. Troubleshooting Guide

Problem	Cause	Solution
The Power Button nor Power Bank LED light up when: - the Power Bank is connected to the heating belt and - the Power Button is pressed.	The Power Bank has been depleted.	Charge the Power Bank.
	The Safety Shut-off Sensor was triggered because of overheating.	Allow the heating belt to cool down to room temperature and switch on again.

The Power Button nor Power Bank LED light up when: - the fully charged Power Bank is connected to the heating belt and - the Power Button has been pressed.	The heating belt or the Power Bank is faulty.	Contact Customer Service for assistance.

8. Technical data

Power Bank:	Type WE-213
Capacity	3.6V ⎓2.5Ah battery capacity, 9Wh
Battery type	Lithium ion 18650 battery
Input / Output Voltage	5V/1A
Operating temperature	14 °F to 104 °F (-10 °C to +40 °C)
Storage temperature	-4 °F to 113 °F (–20 °C to +45 °C)
Recommended storage temperature	14 °F to 77 °F (-10 °C to +25 °C)
Product size	3.7 x 1.2 x 0.8 in (93 x 30 x 21 mm)
Weight	2.4 oz (67 g

For technical data, see the type plate label on the heating belt and on the Power Bank.
USB® is a registered trademark of USB Implementers Forum Inc.

9. FCC Compliance Statement

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

— Reorient or relocate the receiving antenna.
— Increase the separation between the equipment and receiver.
— Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
— Consult the dealer or an experienced radio/TV technician for help.

10. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer Mobile Heating Belt To Go Model HK 67, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.
We will, at our option, repair or replace the Beurer Mobile Heating Belt To Go, Model HK 67, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.
For warranty service contact our customer service department at **1-800-536-0366 or at info@beurer.com** to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase,

your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.
This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment.
The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.
This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.
THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.
IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.
Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.
This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.
This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

Distributed by:
Beurer North America, LP
900 N Federal Hwy, Ste 300
Hallandale Beach, FL 33009 USA

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USB® est une marque commerciale déposée de USB Implementers Forum Inc.