All-in-One Total Solution
Fingers Up for Health !

3in1



EasyTouch®


Model:ET-321

Blood Glucose/Cholesterol/ Hemoglobin Monitoring System


Bioptik is a pillar of innovation in Taiwan, a publicly listed company respected worldwide for its advancements in medical devices.

Step Instructions

Step1.
Check code number and insert code key.

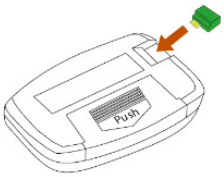

Step2.
Insert test strip and confirm code number.

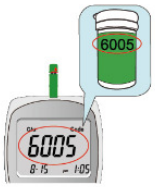

Step3.
Collect blood sample.



⚠ Depends on different types of test strips, the blood might absorb from the top or the side.

3 IN 1	Strip Expiration Detection	Small Blood Sample Needed	Quick Results	Automatic Identification
C G Hb		4 5cm		AUTO

Function		Glucose	Hemoglobin	Cholesterol
Strip Model		SG119 	SH102 	SC142 
Measuring Range		20~600mg/dL (1.1~33.3mmol/L)	7~26 g/dL (4.3~16.1 mmol/L)	100~400mg/dL (2.6~10.4mmol/L)
Memory Capacity		200 results	50 results	50 results
Test Time		10 Sec.	6 Sec.	15 Sec.
Operating Conditions		14~40°C(57~ 104 °F)，20%~85% RH		
Storage Conditions	Meter	-10 ~ +60 °C (14 ~ 140 °F) ，≤95% R.H.		
	Strip	4 ~ 30 °C (39 ~ 86 °F)，≤95% R.H.		
Sample Volume		≧ 4.0 µl	≧ 2.6 µl	≧ 3.0 µl
Sample Type		Finger capillary whole blood		
Strip Shelf Life		24 month	24 month	24 month
Precision		SD≦7.5mg/dL (≦100mg/dL) CV≦6% (>100mg/dL)	CV≦8% (≦9 g/dL) CV≦6% (>9 g/dL)	CV≦8% (≦140mg/dL) CV≦6% (>140mg/dL)
Accuracy		95% within ±15mg/dL (≦100 mg/dL) 95% within ±15% (>100 mg/dL)	95% within ±20%	
Linearity Study		>0.99	>0.9	>0.9
Hematocrit Range		35~50%	20~70%	30~50%
Battery		1.5V (AAA)*2		
Dimensions		88*64*22 H*W*D (mm)		
Display		LCD (35*45 mm)		
Weight		59 g without battery		

Website

Quality Certification

ISO 13485 : 2016 CERTIFIED

QMS

CE 0197

Bioptik EMPOWERING YOU TO LIVE WELL

BIOPTIK TECHNOLOGY, INC. Headquarters
No. 188, Jhonghua S. Rd., Gongguan Village,
Jhunan Township, Miaoli County, 35057, Taiwan