EasyTouch® GCU

EasyTouch® GCU Blood Glucose / Cholesterol / Uric Acid Multi-Function Monitoring System



Instructions for Puncturer



Unscrew Puncturer's Tip



Insert Lancet



Unscrew Protective Cover



Screw Puncturer's Tip




Adjust The Best Depth of Skin Penetration



Pull Up The Sliding Barrel



Hear "Click" Sound Release The Sleeve

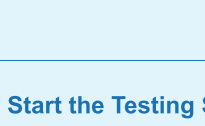

Start the Testing Steps

How to eject the lancet



Unscrew Puncturer's Tip



Push the lancet ejector forward with the thumb and pull out the sliding barrel.

Glucose Testing Step		Cholesterol Testing Step	Uric Acid Testing Step
	Check Code Number		
	1 Insert Code Card		
	Turn On Automatically Insert Test Strip		
	2  Collect Blood Sample		Wipe away the first drop of blood in you finger, and apply the second one to the test strip.
	3 Absorption	 Platform	
	4 Display Result	 Platform	
6 seconds		150 seconds	6 seconds

SETTING TIME & DATE

The meter will enter setting mode automatically when you insert batteries. Please set correct time and date before you begin testing.

Slide battery cover off and insert two alkaline batteries(AAA, 1.5Vx2). into the battery compartment.



1



The meter will turn on automatically.



2



The month will flash.
Press『S』button to set the correct month.



3



Press『M』button to move to the date setting.
The date will flash.
Press『S』button to set the correct date.

4

Press『M』button to move to the hour.
The hour will flash.
Press『S』button to set the correct hour.



5



Press『M』button to move to the minuts.
The minute will flash.
Press『S』button to set the correct minute.



6



The meter will automatically turn off after 30 seconds, while “minute” is set.

CE 0197

Rev.A 12/10
P/N 202336