IM1-HEM-7183-CAP/HEM-7183-LAP-01-06/2024
3150319-6B

OMRON

Automatic Blood Pressure Monitor

Model EZ Power (HEM-7183-C)
EZ Large (HEM-7183-L)


English

Instruction Manual 1: Safety and Other Information

For symbols information, refer to "Symbols Description" section.

Intelli sense

CE 0197

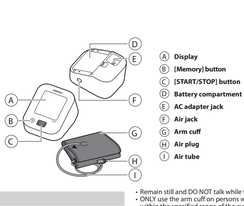Introduction

Thank you for purchasing the OMRON Automatic Blood Pressure Monitor. This blood pressure monitor uses the oscillometric method of blood pressure measurement. This means this monitor detects your blood movement through your brachial artery and converts the movements into a digital reading.

Safety Instructions

This instruction manual provides you with important information about the OMRON Automatic Blood Pressure Monitor. To ensure the safe and proper use of this monitor, READ and UNDERSTAND all of these instructions. **If you do not understand these instructions or have any questions, contact your local OMRON representative before attempting to use this monitor. For specific information about your own blood pressure, consult with your physician.**

Intended Use

Intended Purpose

This device is a digital monitor intended for use in measuring blood pressure and pulse rate in adult patient population. The device also detects the appearance of irregular heartbeats during measurement and indicates this via a symbol with the measurement result.

Intended Patient

Adult patient population

Intended Users

Adult population (may include patients themselves) who can understand this instruction manual.

Clinical Benefit

Patient's blood pressure can be measured non-invasively and simply in the home environment.

Type of Use

This monitor is intended to be multiple patient multiple use.

Limitation

Patient's arm circumference must be 17 - 42 cm.

Indication

This device is used by healthy individuals, patients with hypertension, patients with health-conscious individuals, in a general household situation for measuring blood pressure and pulse rate.

Receiving and Inspection

Remove this monitor and other components from the packaging and inspect for damage. If this monitor or any other component is damaged, DO NOT USE and consult with your local OMRON representative.

Important Safety Information

Read the Important Safety Information in this instruction manual before using this monitor. Follow this instruction manual carefully for your safety.

Keep for future reference. **For specific information about your own blood pressure, consult with your physician.**

Contraindications

- DO NOT use this monitor on an injured arm or an arm under medical treatment.
- DO NOT apply the arm cuff on your arm while on an intravenous drip or blood transfusion.
- DO NOT use this monitor on infants, toddlers, children or persons who cannot express themselves.

Side Effects

- Taking measurements more often than necessary may cause bruising due to blood flow interference.
- Inflating to a higher pressure than necessary may result in bruising of the arm where the cuff is applied. NOTE: refer to "If your systolic pressure is more than 210 mmHg" in the end of instruction manual 2 for additional information.
- Stop using this monitor and consult with your physician if you experience skin irritation or discomfort.

Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.

- DO NOT adjust medication (including altering your dose of any drug or treatment) based on readings from this blood pressure monitor. Take medication as prescribed by your physician. Only a physician is qualified to diagnose and treat high blood pressure and heart related conditions.
- NEVER diagnose or treat yourself based on your readings. ALWAYS consult with your physician.
- If you are experiencing any symptoms or concerns, consult with your physician.
- DO NOT poitpone/stop regular checkups or physician visits based on the results that you are getting from this monitor.
- DO NOT use this monitor in areas containing high frequency (HF) surgical equipment, magnetic resonance imaging (MRI) equipment, computerized tomography (CT) scanners. This may result in incorrect operation of the monitor and/or cause an inaccurate reading.
- DO NOT use this monitor in oxygen rich environments or near flammable gas.
- Consult with your physician before using this monitor if you have common arrhythmias such as atrial or ventricular premature beats or atrial fibrillation; arterial sclerosis; poor perfusion; diabetes; pregnancy; pre-eclampsia or renal disease. NOTE that any of these conditions in addition to patient motion, tremors or shivering may affect the measurement reading.
- To help avoid strangulation, keep the air tube and AC adapter (optional accessory for EZ Large (HEM-7183-L)) cable away from infants, toddlers and children.
- This product contains small parts that may cause a choking hazard if swallowed by infants, toddlers and children.

AC Adapter (optional accessory for EZ Large (HEM-7183-L)) Handling and Usage

- DO NOT use the AC adapter if this monitor or the AC adapter cable is damaged. If this monitor or the cable is damaged, turn off the power and unplug the AC adapter immediately.
- Plug the AC adapter into the appropriate voltage outlet. DO NOT use in a multi-outlet plug.
- NEVER plug in or unplug the AC adapter from the electric outlet with wet hands.
- DO NOT disassemble or attempt to repair the AC adapter.

Battery Handling and Usage

- Keep the batteries out of the reach of infants, toddlers and children.

Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury to the user or patient, or cause damage to the equipment or other property.

- Consult with your physician before using this monitor on an arm where intravascular access or therapy, or an arterio-venous (A-V) shunt, is present because of temporary interference to blood flow and could result in injury.
- Consult with your physician before using this monitor if you have had a mastectomy or lymph node clearance.
- Consult with your physician before using this monitor if you have severe blood flow problems or blood disorders as cuff inflation can cause bruising.
- ONLY inflate the arm cuff when it is applied on your upper arm.
- Remove the arm cuff if it does not start deflating during a measurement.
- DO NOT use this monitor for any purpose other than measuring blood pressure.
- During measurement, make sure that no mobile device or any other electrical device that emits electromagnetic fields is within 30 cm of this monitor. This may result in incorrect operation of the monitor and/or cause an inaccurate reading.
- DO NOT disassemble or attempt to repair this monitor or other components. This may cause an inaccurate reading.
- DO NOT use in a location where there is moisture or a risk of water splashing this monitor. This may damage this monitor.
- DO NOT use this monitor in a moving vehicle such as in a car or on an aircraft.
- DO NOT drop or subject this monitor to strong shocks or vibrations.
- DO NOT use this monitor in places with high or low humidity or high or low temperatures. Refer to section 3.
- During measurement, observe the arm to ensure that the monitor is not causing prolonged impairment to blood circulation.
- DO NOT use this monitor with other medical electrical (ME) equipment simultaneously. This may result in incorrect operation of the devices and/or cause an inaccurate reading.
- Avoid bathing, drinking alcohol or caffeine, smoking, exercising and eating for at least 30 minutes before taking a measurement.
- Rest for at least 5 minutes before taking a measurement.
- Remove tight-fitting or thick clothing from your arm while taking a measurement.
- Remain still and DO NOT talk while taking a measurement.
- ONLY use this arm cuff on persons whose arm circumference is within the specified range of the cuff.
- Ensure that this monitor has acclimated to room temperature before taking a measurement. Taking a measurement after an extreme temperature change could lead to an inaccurate reading. It is recommended that you wait for approximately 2 hours for the monitor to warm up or cool down when the monitor is used in an environment within the temperature specified as operating conditions after it is stored either at the maximum or at the minimum storage temperature. For additional information on operating and storage / transport temperature, refer to section 3.
- DO NOT use this monitor after the durable period has ended. Refer to section 3.
- DO NOT crease the arm cuff or the air tube excessively.
- DO NOT fold or kink the air tube while taking a measurement. This may cause an injury by interrupting blood flow.
- To unplug the air plug, pull on the plastic air plug at the base of the tube, not the tube itself.
- ONLY use the AC adapter (optional accessory for EZ Large (HEM-7183-L)), arm cuff, batteries and accessories specified for this monitor. Use of unsupported AC adapters, arm cuffs and batteries may damage your monitor. See "Optional Accessories".
- ONLY use the approved arm cuff for this monitor. Use of other arm cuffs may result in incorrect readings.
- Read and follow the "Correct Disposal of This Product" in section 4 when disposing of the device and any used accessories or optional parts.

AC Adapter (optional accessory for EZ Large (HEM-7183-L)) Handling and Usage

- Fully insert the AC adapter into the outlet.
- When unplugging the AC adapter from the outlet, be sure to safely pull from the AC adapter. DO NOT pull from the AC adapter cable.
- When handling the AC adapter cable:
 - DO NOT damage it. / DO NOT break it. / DO NOT tamper with it. / DO NOT pinch it. / DO NOT bend it excessively. / DO NOT twist it. / DO NOT use it if it is gathered in a bundle. / DO NOT place it under heavy objects.
- Wipe any dust off of the AC adapter.
- Unplug the AC adapter when not in use.
- Unplug the AC adapter before wiping this monitor.

Battery Handling and Usage

- DO NOT insert batteries with their polarities incorrectly aligned.
- ONLY use 4"AA" alkaline or manganese batteries with this monitor.
- DO NOT use other types of batteries. DO NOT use new and used batteries together. DO NOT use different brands of batteries together.
- Remove the batteries if this monitor will not be used for a long period of time.
- If battery fluid should get in your eyes, immediately rinse with plenty of clean water. Consult with your physician immediately.
- If battery fluid should get on your skin or clothes, wash your skin immediately with plenty of clean, lukewarm water. If irritation, injury or pain persists, consult with your physician.
- DO NOT use batteries after their expiration date.
- Periodically check the batteries to ensure they are in good working condition.

General Notices

- To stop a measurement, press the [START/STOP] button while taking a measurement.
- When you take a measurement on the right arm, the air tube should be at the side of your elbow. Be careful not to rest your arm on the air tube.

- Blood pressure may differ between the right and left arm, and may result in a different measurement value. Always use the same arm for measurements. If the values between both arms differ substantially, check with your physician which arm to use for your measurements.
- When using an AC adapter, make sure not to place your monitor in a location where it is difficult to plug and unplug the AC adapter.

Battery Handling and Usage

- Disposal of used batteries should be carried out in accordance with local regulations.
- The supplied batteries may have a shorter life span than new batteries.
- Replacing batteries will not delete previous readings.

1 Error Messages and Troubleshooting

If any of the below problems occur during measurement, check to make sure that the other electrical devices within 30 cm of the monitor. If the problem persists, refer to the table below.

Display/Problem	Possible Cause	Solution
E1 appears or the arm cuff does not inflate.	The [START/STOP] button was pressed while the arm cuff is not applied.	Press the [START/STOP] button again to turn the monitor off.
	Air plug is not completely inserted.	Insert the air plug securely.
	The arm cuff is not applied correctly.	Apply the arm cuff correctly, then take another measurement. Refer to section 4 of instruction manual 2.
	Air is leaking from the arm cuff.	Replace the arm cuff with a new one. Refer to section 10 of instruction manual 2.
E2 appears or a measurement cannot be completed after the arm cuff inflates.	Moving or talking during a measurement causes the arm cuff to not inflate sufficiently.	Remain still and do not talk during a measurement. If "E2" appears repeatedly, inflate the arm cuff manually until the systolic pressure is 30 to 40 mmHg above your previous readings. Refer to the end of instruction manual 2.
	The systolic pressure is above 210 mmHg and a measurement cannot be taken.	Do not touch the arm cuff and/or bend the air tube while taking a measurement. If inflating the arm cuff manually, refer to the end of instruction manual 2.
E3 appears	The arm cuff is inflated exceeding the maximum allowable pressure.	Moving or talking during a measurement results in vibrations that disrupt the measurement.
E4 appears	Moving or talking during a measurement results in vibrations that disrupt the measurement.	Remain still and do not talk during a measurement.
E5 appears	The pulse rate is not detected correctly.	Apply the arm cuff correctly, then take another measurement. Refer to section 4 of instruction manual 2. Remain still and sit correctly during a measurement.
Er appears	The monitor has malfunctioned.	Press the [START/STOP] button again. If "Er" still appears, contact your OMRON retail outlet or distributor.
Display flashes	The pulse rate is not detected correctly.	Apply the arm cuff correctly, then take another measurement. Refer to section 4 of instruction manual 2. Remain still and sit correctly during a measurement.
Batteries are low.	Replacing all 4 batteries with new ones is recommended. Refer to section 3 of instruction manual 2.	
and appears, or the monitor is turned off unexpectedly during a measurement.	Batteries are depleted.	Immediately replace all 4 batteries with new ones. Refer to section 3 of instruction manual 2.
Nothing appears on the display of the monitor.	Battery polarities are not properly aligned.	Check the battery installation for proper alignment. Refer to section 3 of instruction manual 2.
Readings appear too high or too low.	Blood pressure varies constantly. Many factors including stress, time of day, and/or how you take the air, may affect your blood pressure. Review sections 2, 4 and 5 of instruction manual 2.	
Any other problem occurs.	Press the [START/STOP] button to turn the monitor off, then press it again to take a measurement. If the problem continues, remove all batteries and wait for 30 seconds. Then re-install the batteries. If the problem still persists, contact your OMRON retail outlet or distributor.	

2 Maintenance

2.1 Maintenance

To protect your monitor from damage, follow the directions below: Changes or modifications not approved by the manufacturer will void the user warranty.

⚠ Precaution

DO NOT disassemble or attempt to repair this monitor or other components. This may cause an inaccurate reading.

2.2 Storage

- Store your monitor and other components in a clean, safe location.
- Gently fold the air tube into the arm cuff. Note: Do not bend or crease the air tube excessively.
- Do not store your monitor and other components:
 - If your monitor and other components are wet.
 - In locations exposed to extreme temperatures, humidity, direct sunlight, dust or corrosive vapors such as bleach.
 - In locations exposed to vibrations or shocks.

2.3 Wiping the Monitor

- Do not use any abrasive or volatile cleaners.
- Use a soft dry cloth or a soft cloth moistened with mild (neutral) detergent to wipe your monitor and arm cuff, and then wipe them with a dry cloth.
- Do not wash or immerse your monitor and arm cuff or other components in water.
- Do not use gasoline, thinners or similar solvents to wipe your monitor and arm cuff or other components.

2.4 Calibration and Service

- The accuracy of this blood pressure monitor has been carefully tested and is designed for a long service life.
- It is generally recommended that the unit inspected every two years to ensure correct functioning and accuracy. Please consult your authorised OMRON dealer or the OMRON Customer Service at the address given on the packaging or attached literature.

3 Specifications

Product description	Automatic Blood Pressure Monitor
Model	EZ Power (HEM-7183-C)/ EZ Large (HEM-7183-L)
Display	LCD digital display
Cuff pressure range	0 to 299 mmHg
Blood pressure measurement range	SYS: 60 to 260 mmHg DIA: 40 to 215 mmHg
Pulse measurement range	40 to 180 beats / min.
Accuracy	Pressure: ±3 mmHg Pulse: ±5 % of display reading
Measurement method	Oscillometric method
Operation mode	Continuous operation
IP classification	Monitor: IP21 AC adapter: IP21 (HHP-CM01)/IP22 (HHP-CM01, HHP-CH01)
Rating	DC6 V 4 W
Power source	4 "AA" batteries 1.5 V or AC adapter (INPUT: AC 100 - 240 V 50 - 60 Hz 0.12 - 0.065 A)
Battery life	Approx. 1500 measurements (Using new alkaline batteries and included arm cuff. Depending on the type of battery and conditions.)
Durable period (Service life)	Monitor: 5 years or the time when reaching 30000 times of use. / Cuff: 5 years or the time when reaching 10000 times of use. / AC adapter: 5 years
Operating conditions	+10 to +40 °C / 15 to 90 % RH (non-condensing) / 800 to 1060 hPa
Storage / Transport conditions	-20 to +60 °C / 10 to 90 % RH (non-condensing) / 500 to 1060 hPa
Weight	Monitor: approx. 260 g (not including batteries) Arm cuff for EZ Power (HEM-7183-C): approx. 110 g Arm cuff for EZ Large (HEM-7183-L): approx. 170 g
Dimensions	Monitor: approx. 98 mm (W) × 79 mm (H) × 138 mm (L) / Arm cuff for EZ Power (HEM-7183-C): approx. 146 mm × 466 mm (air tube: 610 mm) / Arm cuff for EZ Large (HEM-7183-L): approx. 145 mm × 594 mm (air tube: 750 mm)
Cuff circumference applicable to the monitor	17 to 42 cm (included arm cuff for EZ Power (HEM-7183-C): 22 to 32 cm / for EZ Large (HEM-7183-L): 22 to 42 cm)
Memory	Stores up to 30 readings
Contents	Monitor, arm cuff (HEM-CR24 for EZ Power (HEM-7183-C) / HEM-RML31 for EZ Large (HEM-7183-L)), AC adapter (HHP-CM01 or HHP-CH01: included only in EZ Power (HEM-7183-C)), "4 A&D batteries", Instruction Manual 1 and 2, warranty card
Protection against electric shock	Internally powered ME equipment (when using only batteries) Class II ME equipment (when using AC adapter)
Applied part	Type BF (arm cuff)

Note

- These specifications are subject to change without notice.
- This monitor is investigated according to the clinical requirements of EN ISO 81060-2:2014 and complies with EN ISO 81060-2:2019+A1:2020+A2:2024. In the clinical validation study, KS was used on 85 subjects for determination of diastolic blood pressure.
- This monitor has not been validated for use on pregnant patients.
- IP classification is degrees of protection provided by enclosures in accordance with IEC 60529. This monitor and AC adapter (included only in EZ Power (HEM-7183-C)) are protected against solid foreign objects of 12.5 mm diameter and greater such as a finger. The monitor and AC adapter HHP-CM01 (included only in EZ Power (HEM-7183-C)) are protected against vertically falling water drops which may cause issues during a normal operation. The AC adapter HHP-CH01 and HHP-CM01 are protected against oblique falling water drops which may cause issues during a normal operation.
- Operation mode classification complies with IEC 60601-1.
- Please refer to "UDI [8012].x,x" for the software version/build number; "x,x,x" may consist of a set of 2 to 6 numbers.

4 Correct Disposal of This Product (Waste Electrical & Electronic Equipment)

This marking shown on the product or its literature, indicates that it should not be disposed of, with other household waste at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this product from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take this item for environmentally safe recycling.
Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial waste for disposal.

5 Important Information regarding Electromagnetic Compatibility (EMC)

This device conforms to the EN 60601-1-2:2015+A1:2021 Electromagnetic Compatibility (EMC) standard.
Further documentation in accordance with this EMC standard is available at:
https://healthcare.omron.com/electro-magnetic-compatibility
Special precautions for this device on the website.

6 Guidance and Manufacturer's Declaration

- This OMRON product is produced under the strict quality system of OMRON HEALTHCARE Co., Ltd., Japan. The Core component for OMRON blood pressure monitors, which is the Pressure Sensor, is produced in Japan.
- Please report to the manufacturer and the regulatory authority in which you are established about any serious incident that has occurred in relation to this device.

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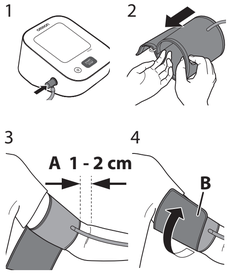
A. Tube side of the cuff should be 1 - 2 cm above the inside of your elbow.
B. Make sure that the air tube is on the inside of your arm and wrap the cuff securely so it can no longer slip round.

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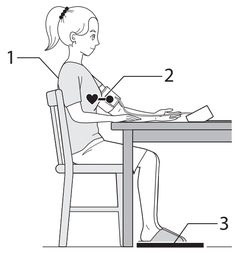
Relax and sit comfortably. Remain still and do not talk.
1. Keep your back and arm supported.
2. Keep the arm cuff at the same level as your heart.
3. Keep your feet flat and your legs uncrossed.

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
6 Taking a Measurement

TH การทำการวัด
ID Melakukan Pengukuran
VI Tiến hành đo
MS Mengambil Ukuran

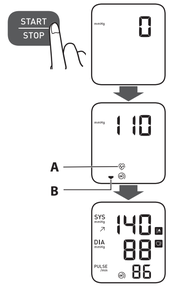1. Press the [START/STOP] button.
A: Flashes at every heartbeat.
B: Appears while the cuff is deflating.
2. The reading is saved automatically.

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E1
E2
• • •

Instruction Manual 1
2.

Instruction Manual 2: Operational Instructions

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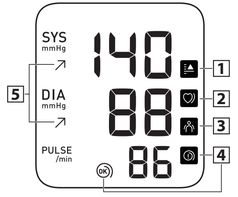
*** The high blood pressure definition is based on the 2023 ESH Guideline.**
**** An irregular heartbeat rhythm is defined as a rhythm that is 25 % less or 25 % more than the average rhythm detected while your monitor is measuring blood pressure.**

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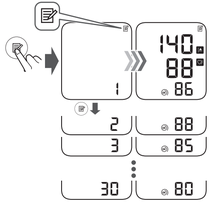
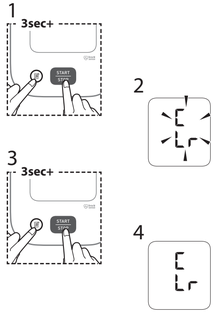
Model: HEM-CS24
17 - 22 cm

AC Adapter

Model: HHP-CM01
HHP-OH01
(For Australia)
(Only for EZ Large (HEM-7183-L))

Manufacturer	ผู้ผลิต Produsen	Nhà sản xuất Pengilang	**OMRON HEALTHCARE Co., Ltd.** 53, Kunotsubo, Terado-cho, Muko, KYOTO, 617-0002 JAPAN
EC REP EU-representative	ตัวแทนที่ได้รับมอบอำนาจในยุโรป (EU-representative) Perwakilan UE	Đại diện - EU Wakil EU	**OMRON HEALTHCARE EUROPE B.V.** Wegalaan 73, 2132 JD Hoofddorp, THE NETHERLANDS www.omron-healthcare.com
Importer in EU	ผู้นำเข้าในสหภาพยุโรป Importir di UE	Nhà nhập khẩu tại EU Pengimport di EU	
Asia Pacific HQ	สำนักงานใหญ่ในเอเชียแปซิฟิก Kantor Pusat Asia Pasifik	Trụ sở Khu vực Châu Á Thái Bình Dương Ibu Pejabat Asia Pasifik	**OMRON HEALTHCARE SINGAPORE PTE LTD.** 438B Alexandra Road, #08-01/02 Alexandra Technopark, Singapore 119968 www.omronhealthcare-ap.com
Production Facility	โรงงานผลิต Fasilitas Produksi	Cơ sở sản xuất Kemudahan Pengeluaran	**OMRON HEALTHCARE MANUFACTURING VIETNAM CO., LTD.** No.28 VSIP II, Street 2, Vietnam-Singapore Industrial Park II, Binh Duong Industry-Service-Urban Complex, Hoa Phu Ward, Thu Dau Mot City, Binh Duong Province, Vietnam

Issue Date / วันที่ออกจำหน่ายฉบับ / Tanggal Terbit / Ngày phát hành / Tarikh Keluaran : 2024-12-02

Made in Vietnam /ผลิตในเวียดนาม / Buatan Vietnam / Sản xuất tại Việt Nam / Buatan Vietnam

https://www.omronhealthcare-ap.com/