GEA medical

OPERATOR'S MANUAL



YC-B INFANT BED

DEPKES RI AKL 2090 2901 182

Ⅰ. APPLICATION


Infant bed is suitable for observing, suckling or nursing the neonates in infant room, delivery room and department of paediatrics.


Ⅱ. COMPOSITION OF PRODUCTS


The infant bed consists of bed, bed frame, mattress, and I.V. pole.

Ⅲ. ASSEMBLING

Generally, the Infant bed is packed to one carton with the figuration size of 910mm×420mm×800mm. Take care not to scratch or otherwise damage the surfaces of the frame when unpacking. As following figures indicate, install the wheel and I.V. pole.

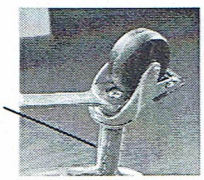


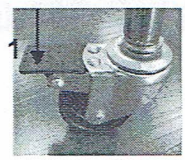


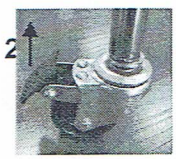

Put on the spring washer 10 and flat washer 10 on the wheels, and then connect them to bed frame with special spanner firmly.

Note: lock the wheels as the arrow 1 indicates, and unlock the wheels as the arrow 2 indicates.

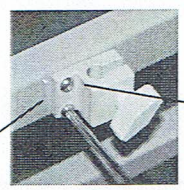

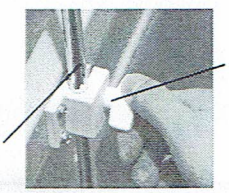

Insert cross groove bolt M5×20 through Seat for I.V. Pole and bed frame, then tighten them with flat washer 5, spring Washer 5 and hexagon nut M5.

Insert the I.V. pole into the fixed seat, then adjust it to a proper height and screw down the tightening screw.

Ⅳ. SYMBOLS

Date Of Manufacture

SN Serial Number

Ⅴ. SPECIFICATION

1. The Max. Load of bassinet is 5 Kg.
2. The Max. Load of I.V. pole is 2 Kg.
3. Toxicity for the material which contact with the patient's skin directly is less than Class I. No skin stimulus and delayed hypersensitivity as well.

Ⅵ. OPERATION

The Infant bed can be just placed by side of mother's bed as the wheels must be locked. Unlock the wheels of the infant bed when it needs to be moved. A drawer attached to the bed frame is used to store some baby's clothing or diapers and so on.

Ⅶ. NOTICE

1. Don't overload the bassinet, or else it will be damaged.
2. Please ensure the tightening screw is firmly screwed down, or else the I.V. pole might fall off.
3. Do not use alcohol, acetone, or any organic solvents for cleaning the plexiglass. Do not expose the plexiglass to direct ultraviolet radiation as well.

Ⅷ. THE LIFETIME OF PRODUCT

Its lifetime is 3 years, lifetime means the period from sell-by date to the date of discarding as useless. When the lifetime of the product ends, please stop using it and handle it properly.

Ⅸ. CLEANING AND MAINTENANCE

When first operate the head frame or the operation to infant with the head frame over, the head frame should be cleaned and disinfected drastically. When the operation of head frame is over each time, it should be cleaned and disinfected too. Please use a disinfectant-detergent (like 84 detergent) to thoroughly clean all surfaces, including all holes and indentations. Unpack the bed cover of mattress and use detergent for cleaning. Dry the bed cover after it cleaned in the water.

The infant bed should be stored in dry and cool place with the relative humidity ≤80% and without any corrosive gas after operation.

Ⅹ. PACKING LIST

NO.	DESCRIPTION OF PART	QUANTITY	NO.	DESCRIPTION OF PART	QUANTITY
1	Bed	1	2	Bed frame	1
3	Mattress	1	4	I.V. pole (including fixed seat and tightening screw)	1
5	Wheel with brake	4	6	Special spanner	1
7	Certificate	1			