Sinocare

BSX-516



USB Port Support

WHO Function

Large LCD Display

Back-Light Function

Automatic Turn-off

One-year Warranty

Oscillometric Method


Two Users Memory

Add:No.265 Guyuan Road Hi-tech Zone,Changsha,Hunan China,410205
Tel:0731-89935929 Email:intl_mkt@sinocare.com Web:www.sinocareintl.com


Technical Specifications

Product Type:	BSX516
Cuff Circumference:	22cm~42cm(±5mm)
Memory:	2*99 groups memory (Double users)
Measurement Range:	Blood pressure: 0-290 mmHg Pulse 40-180 beats / min
Accuracy:	Blood pressure : + 3 mmHg Pulse : ±5% of the reading
Material:	ABS + LCD digital display
Power Source:	Internal-DC 6V (4 * AA) / External-DC6V 0.6A
Operating Conditions:	Temperature: 5 °C to 40 °C Humidity: 15% to 85% HR Pressure altitude: 80Kpa ~ 106Kpa
Storage And Shipping Conditions:	Temperature: -20 °C a 55 °C Humidity: ⩽93% RH Pressure altitude: 70Kpa ~ 106Kpa
Measurement Mode:	Oscillometric method


Use methods and precautions

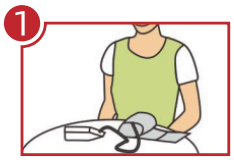


Keep relaxed before measuring sit down quietly for a moment

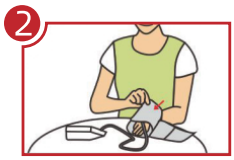

Palms up, put the arm cuff parallel to the heart

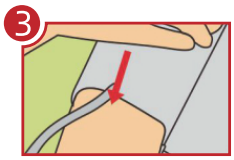


Palms up, keep the intake pipe and the arteries parallelled

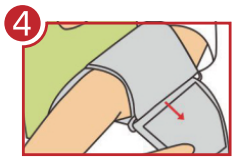


Wrap the arm cuff around your arm tightly in the opposite direction, paste together



Place one finger into it is the most appropiate

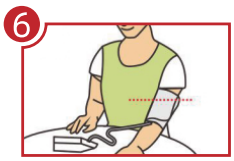

Keep the arm cuff parallelled to the heart, palms up

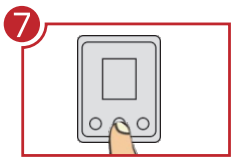

Press the on/off button, keep reaxed and start measuring



High pressure, low pressure, pulse instantly displayed after 40 seconds